



## **Dioxins and benzofurans**

# Analysis of dioxins and benzofurans via comprehensive GC x GC

## Application Note

Environmental

### **Authors**

Agilent Technologies, Inc.

### **Introduction**

Analysis of dioxins and benzofurans by GC, with a PDMS phase equivalent to Agilent VF-1ms, and Agilent VF-23ms, is accomplished in 90 minutes.




**Agilent Technologies**

## Conditions

Technique : GC x GC  
Column 1 : 100% methyl PDMS, Agilent equivalent VF-1 ms, 0.25 mm x 30 m (df = 0.25  $\mu$ m) (Part no. CP8912)  
Column 2 : Agilent VF-23ms, 0.10 mm x 1 m (df= 0.1  $\mu$ m) (Part no. CP8819) as 10 m column  
Temperature 1 : 90 °C, (2 min), 30 °C/min  $\rightarrow$  170 °C, 5 min, 1.5 °C/min  $\rightarrow$  290 °C  
Temperature 2 : 90 °C, (2 min), 30 °C/min  $\rightarrow$  170 °C, 5 min, 1.5 °C/min  $\rightarrow$  290 °C  
Carrier Gas : Helium  
Pressure program : 300 kPa  
Injector : Splitless, T = 290 °C  
Detector : ECD, T = 300 °C  
Odulator : single jet  
Sample Size : 1  $\mu$ L  
Sample : 100 pg/ $\mu$ L

Courtesy : Peter Korytar, The Netherlands Institute for Fisheries Research, Ijmuiden

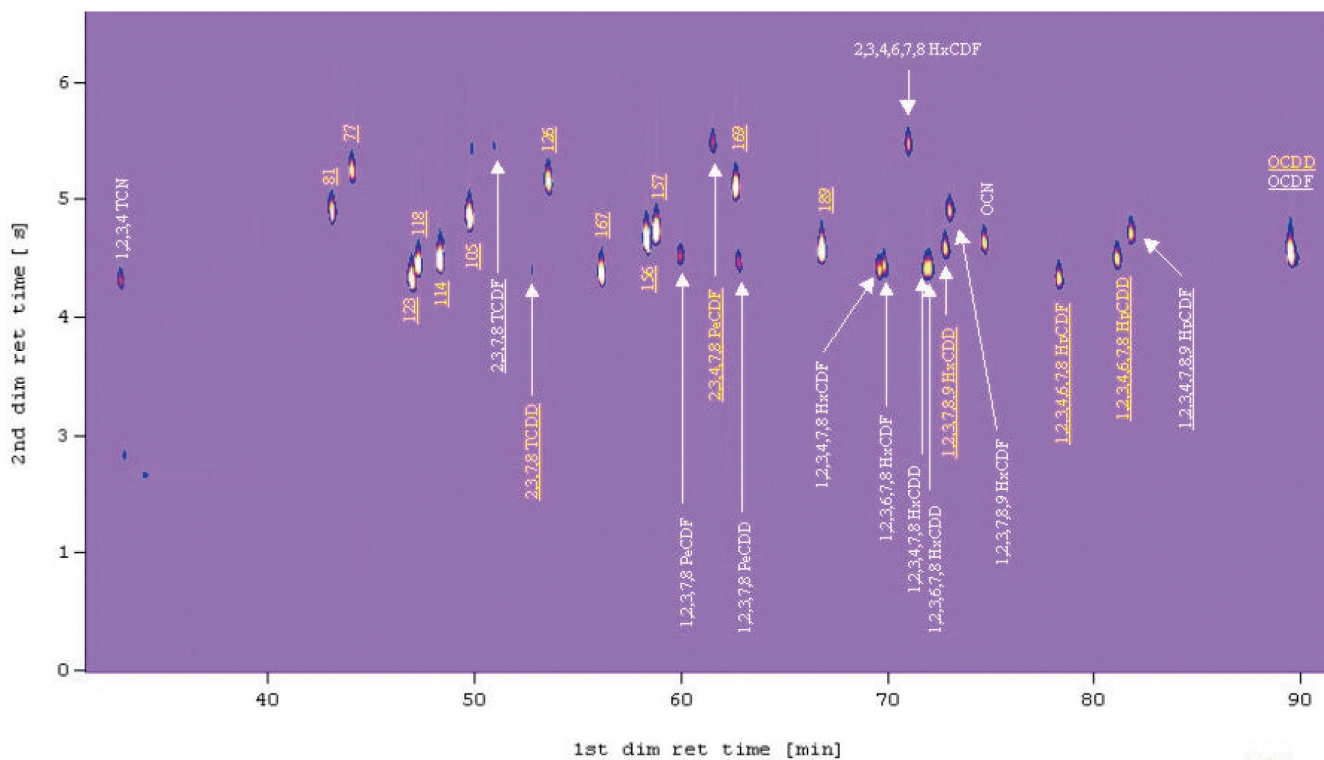

[www.agilent.com/chem](http://www.agilent.com/chem)

This information is subject to change without notice.

© Agilent Technologies, Inc. 2011

Printed in the USA

31 October, 2011

First published prior to 11 May, 2010

A02096



**Agilent Technologies**