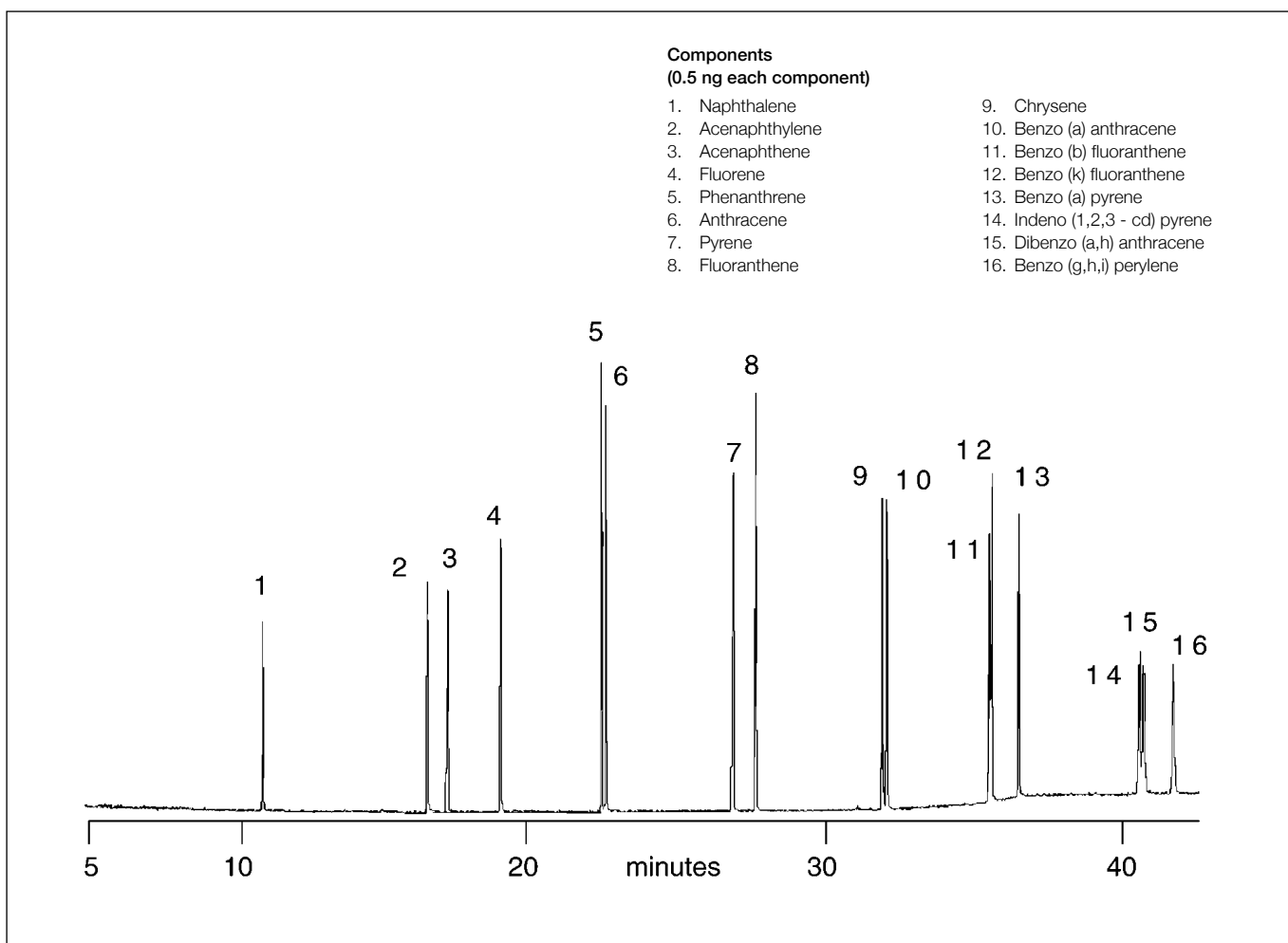


## Analysis of polynuclear aromatic hydrocarbons (PAH)

**BPX5**

Column part number	054113		
Phase	BPX5	Final temperature	290°C, 10 min
Column	25 m x 0.22 mm x 0.25 µm	Detector	MS
Initial temperature	50°C, 2 min	Injection mode	Split (40:1)
Rate	8°C/min	Carrier gas	He, 15 psi


 For more information visit [www.trajanscimed.com](http://www.trajanscimed.com) or contact [techsupport@trajanscimed.com](mailto:techsupport@trajanscimed.com)